

Low-input Capture-C: A Chromosome Conformation Capture Assay to Analyze Chromatin Architecture in Small Numbers of Cells

A. Marieke Oudelaar*, Damien J. Downes, James O.J. Davies and Jim R. Hughes*

Medical Research Council (MRC) Molecular Haematology Unit, Weatherall Institute of Molecular Medicine, University of Oxford, Oxford, United Kingdom

*For correspondence: oudelaar@well.ox.ac.uk; jim.hughes@imm.ox.ac.uk

[Abstract] Chromosome conformation capture (3C) techniques are crucial to understanding tissue-specific regulation of gene expression, but current methods generally require large numbers of cells. This protocol describes two new low-input Capture-C approaches that can generate high-quality 3C interaction profiles from 10,000-20,000 cells, depending on the resolution used for analysis.

Keywords: Nuclear organization, Chromatin architecture, Chromosome conformation capture (3C), Capture-C, Low input, Small cell numbers

[Background] 3C techniques play a key role in investigating how nuclear organization and structural interactions between regulatory elements relate to gene activity (Dekker *et al.*, 2002). As these interactions are highly tissue-specific, it is crucial that 3C experiments are performed in well defined, purified cell populations.

A major limitation of 3C techniques is the large numbers of cells required: current methods use between 100,000 and 10 million cells (Davies *et al.*, 2017). Many primary tissues and rare cell populations are not available in these numbers. We have therefore developed two new low-input Capture-C approaches that can generate high-quality interaction profiles from ~20,000 cells at maximum resolution (individual DpnII fragments), and from ~10,000 cells using windowing based analysis (Oudelaar *et al.*, 2017).

Materials and Reagents

A. 3C library preparation

1. LowBind DNA 1.5 ml microtube
2. Phase Lock Gel (PLG) Light tubes (VWR, Quantabio, catalog number: 733-2477)
3. Genomic DNA ScreenTape (Agilent Technologies, catalog number: 5067-5366)
4. TapeStation Loading Tips (Agilent Technologies, catalog number: 5067-5153)
5. 1 M glycine (Sigma-Aldrich, catalog number: G7126)
6. Phosphate-buffered saline (PBS) (Thermo Fisher Scientific, Gibco™, catalog number: 10010031)
7. Growth media
8. 37% formaldehyde (Sigma-Aldrich, catalog number: 252549)

9. Ethanol absolute > 99.8% (VWR, catalog number: 20821.330, or equivalent)
10. Dry ice or liquid nitrogen
11. PCR grade water (Thermo Fisher Scientific, Invitrogen™, catalog number: AM9932)
12. SDS (20% v/v in water) (Thermo Fisher Scientific, Invitrogen™, catalog number: AM9820)
13. Triton X-100 (20% v/v in water) (Sigma-Aldrich, catalog number: T8787)
14. *DpnII* (50,000 U/ml) (New England Biolabs, catalog number: R0543M)
15. T4 DNA ligase (30 U/μl) (Thermo Fisher Scientific, Thermo Scientific™, catalog number: EL0013)
16. Proteinase K (Thermo Fisher Scientific, Thermo Scientific™, catalog number: EO0491)
17. RNase (DNase-free) (Roche Diagnostics, catalog number: 11119915001)
18. Phenol-chloroform-isoamyl alcohol (PCI) 25:24:1 (Sigma-Aldrich, catalog number: 77617)
19. 3 M NaOAc (Thermo Fisher Scientific, Invitrogen™, catalog number: AM9740)
20. GlycoBlue (Thermo Fisher Scientific, Invitrogen™, catalog number: AM9515)
21. Genomic DNA Reagents (Agilent Technologies, catalog number: 5067-5365)
22. TaqMan Universal PCR Master Mix II without UNG (Thermo Fisher Scientific, Applied Biosystems™, catalog number: 4440040)
23. KAPA Sybr Fast Universal (Kapa Biosystems, catalog number: KK4602)
24. Fresh lysis buffer (see Recipes)
 - a. Tris pH 8 (Thermo Fisher Scientific, Invitrogen™, catalog number: AM9855G)
 - b. Sodium chloride (NaCl) (Thermo Fisher Scientific, Invitrogen™, catalog number: AM9760G)
 - c. Igepal CA-630 (100 μl of 10%) (Sigma-Aldrich, catalog number: I8896)
 - d. cOmplete Protease Inhibitor Cocktail (1 tablet in 2 ml of PCR grade water for 25x; store at -20 °C) (Sigma-Aldrich, Roche Diagnostics, catalog number: 11873580001)
 - e. PCR grade water
25. qPCR primers ([Appendix I](#))

B. Addition of Illumina sequencing adaptors by ligation or tagmentation to generate Capture-C libraries

LI-Capture-C

1. Covaris microTUBE AFA Fiber pre-split snap-cap 6 x 16 mm (Covaris, catalog number: 520045)
2. LowBind DNA 1.5 ml microtube
3. PCR tube
4. D1000 ScreenTape (Agilent Technologies, catalog number: 5067-5582)
5. TapeStation Loading Tips (Agilent Technologies, catalog number: 5067-5153)
6. NEBNext Ultra DNA Library Prep Kit for Illumina (New England Biolabs, catalog number: E7370S/L)
7. Agencourt Ampure XP SPRI Beads (Beckman Coulter, catalog number: A63881)
8. NEBNext Multiplex Oligos for Illumina Primer set 1 (New England Biolabs, catalog number: E7335S/L)

9. NEBNext Multiplex Oligos for Illumina Primer set 2 (New England Biolabs, catalog number: E7500S/L)
10. Herculase II Fusion Polymerase Kit (Agilent Technologies, catalog number: 600677)
11. PCR grade water (Thermo Fisher Scientific, Invitrogen™, catalog number: AM9932)
12. Ethanol absolute > 99.8% (VWR, catalog number: 20821.330, or equivalent)
13. D1000 Reagents (Agilent Technologies, catalog number: 5067-5583)
14. Qubit dsDNA BR Assay kit (Thermo Fisher Scientific, Invitrogen™, catalog number: Q32850)

Tag-Capture-C

1. LowBind DNA 1.5 ml microtube
2. PCR tube
3. High Sensitivity D1000 ScreenTape (Agilent Technologies, catalog number: 5067-5584)
4. D1000 ScreenTape (Agilent Technologies, catalog number: 5067-5582)
5. TapeStation Loading Tips (Agilent Technologies, catalog number: 5067-5153)
6. Nextera DNA Sample Preparation Kit (Illumina, catalog number: FC-121-1030/1031)
7. Zymo DNA Clean & Concentrator-5 Kit (Zymo Research, catalog number: D4013)
8. KAPA HiFi PCR Kit with dNTP (KK2102) (Roche Diagnostics, catalog number: 07958846001)
9. Custom designed index primers
10. Agencourt Ampure XP SPRI Beads (Beckman Coulter, catalog number: A63881)
11. PCR grade water (Thermo Fisher Scientific, Invitrogen™, catalog number: AM9932)
12. Ethanol absolute > 99.8% (VWR, catalog number: 20821.330, or equivalent)
13. High Sensitivity D1000 Reagents (Agilent Technologies, catalog number: 5067-5585)
14. D1000 Reagents (Agilent Technologies, catalog number: 5067-5583)
15. Qubit dsDNA BR Assay kit (Thermo Fisher Scientific, Invitrogen™, catalog number: Q32850)

C. Enrichment for viewpoints of interest by oligonucleotide capture

1. Safeseal Microcentrifuge Tubes (Sorenson BioScience, catalog number: 39640T)
2. D1000 ScreenTape (Agilent Technologies, catalog number: 5067-5582)
3. TapeStation Loading Tips (Agilent Technologies, catalog number: 5067-5153)
4. Nimblegen SeqCap EZ Hybridization and wash kit (Roche Diagnostics, catalog number: 05634261001)
5. Nimblegen SeqCap EZ Accessory kit v2 (Roche Diagnostics, catalog number: 07145594001)
6. 1 µg/µl COT DNA of relevant species (Mouse) (Thermo Fisher Scientific, Invitrogen™, catalog number: 18440016)
7. M-270 Streptavidin Dynabeads (Thermo Fisher Scientific, Invitrogen™, catalog number: 65305)
8. Agencourt Ampure XP SPRI Beads (Beckman Coulter, catalog number: A63881)
9. PCR grade water (Thermo Fisher Scientific, Invitrogen™, catalog number: AM9932)
10. Qubit dsDNA BR Assay Kit (Thermo Fisher Scientific, Invitrogen™, catalog number: Q32850)
11. Qubit dsDNA HS Assay Kit (Thermo Fisher Scientific, Invitrogen™, catalog number: Q32851)

12. D1000 Reagents (Agilent Technologies, catalog number: 5067-5583)
13. KAPA Library Quantification Complete Kit (Universal) (Kapa Biosystems, catalog number: KK4824)
14. LI-Capture-C
Nimblegen SeqCap EZ HE-oligo kit A (Roche Diagnostics, catalog number: 06777287001)
Nimblegen SeqCap EZ HE-oligo kit B (Roche Diagnostics, catalog number: 06777317001)
15. Tag-Capture-C
xGen Custom Blocking Oligos (P5/i5 & P7/i7 Nextera blocker)

Equipment

1. Centrifuge
2. Incubator (any)
3. Eppendorf Thermomixer C (or equivalent shaking incubator) (Eppendorf, model: ThermoMixer[®] C, catalog number: 5382000015)
4. P10, P20, P200 and P1000 pipettes (any)
5. Thermocycler (any)
6. NanoDrop
7. Agilent 2200 TapeStation (Agilent Technologies, model: Agilent 2200 TapeStation System)
Note: Assessments on the Agilent TapeStation can also be performed on the Agilent BioAnalyzer.
8. Sonicator (Covaris, model: S220 Focused-ultrasonicator)
9. DynaMag-2 (Thermo Fisher Scientific, catalog number: 12321D)
10. Vacuum Concentrator (SciQuip, CHRIST series)


Procedure

Overview

Note: Overview of the procedure (Graphics adapted from Davies et al., 2017).

A. 3C library preparation (5 days)

1. Fixation—2 h
2. Digestion—1 d
3. Ligation—1 d
4. De-crosslinking—1 d
5. DNA extraction—1 d
6. Quantification and QC—3 h

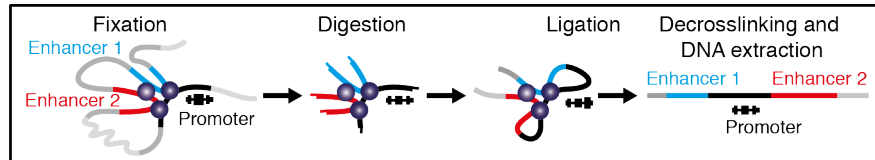

B. Addition of Illumina sequencing adaptors by ligation or tagmentation to generate Capture-C libraries (2 days)


LI-Capture-C:

1. Sonication— $\frac{1}{2}$ d
2. End repair and adaptor ligation— $\frac{1}{2}$ d
3. PCR addition of indices— $\frac{1}{2}$ d

Or


Tag-Capture-C:

1. Tagmentation—1 d
2. PCR addition of indices— $\frac{1}{2}$ d

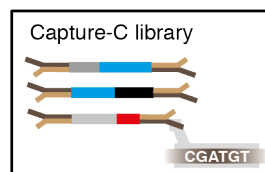

C. Enrichment for viewpoints of interest by oligonucleotide capture (8 days)

1. Capture I
 - a. Hybridization—4 d
 - b. Washing and recovery of captured material—1 d
 - c. PCR amplification of captured material—1 d
2. Capture II
 - a. Hybridization—1 d
 - b. Washing and recovery of captured material— $\frac{1}{2}$ d
 - c. PCR amplification of captured material— $\frac{1}{2}$ d

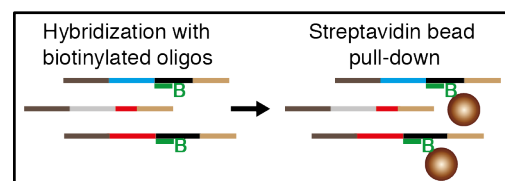

Steps

A. 3C library preparation

1. Fixation (2 h)

Note: The fixation procedure described below has been optimized for erythroid cells. Depending on the cell type used, some of centrifugation steps might need to be adjusted to ensure proper pelleting.

- a. Pre-cool centrifuge to 4 °C. Chill glycine, PBS, and lysis buffer (see Recipes).
- b. Collect cells from tissue and make single-cell suspensions of 10,000 to 10⁶ cells in 1 ml growth media. If sorting cells, sort into 100 µl media and add media afterwards to a total volume of 1 ml.
- c. Add 54 µl 37% formaldehyde (2% final concentration), mix well and incubate for 10 min at RT on a rocking or tumbling incubator.
- d. Quench by adding 150 µl 1 M cold glycine (1/8 of the final volume. *i.e.*, 1.2 ml).
- e. Centrifuge for 15 min/500 x g/4 °C.
- f. Carefully remove supernatant, do not disturb pellet. Can leave ~5% behind.
- g. Wash pellet by gently resuspending in 1 ml cold PBS.
- h. Centrifuge for 15 min/500 x g/4 °C.
- i. Carefully remove supernatant, do not disturb pellet. Can leave ~5% behind.
- j. Resuspend pellet in 500 µl cold lysis buffer.
- k. Incubate for 20 min on ice.
- l. Snap freeze with ethanol and dry ice or liquid nitrogen. Snap freezing aids digestion so cells can be thawed again for digestion at this point or stored long-term at -80 °C.
SAFE STOPPING POINT: Store at -80 °C.

2. Digestion (1 day)

- a. Pre-warm a ThermoMixer to 37 °C.
- b. Prepare Digestion Mix (make a master mix for multiple libraries):

PCR grade water	152.5 µl
10x <i>DpnII</i> Buffer	20 µl
20% SDS	2.5 µl
- c. Defrost aliquots of 10,000 to 10⁶ formaldehyde-fixed cells for each reaction.
- d. Centrifuge for 10 min/15,000 x g/4 °C.
- e. Carefully remove all the lysis buffer, using a P20 to remove residual buffer without disturbing the pellet.
- f. Resuspend in 175 µl of Digestion Mix.
- g. Shake horizontally (1,400 rpm) for 1 h at 37 °C using the ThermoMixer.
- h. Add 16 µl of 20% Triton X-100 (1.7% final concentration) and shake for 1 h at 37 °C.
- i. Add first aliquot of 3 µl of *DpnII* enzyme and shake at 37 °C for at least two hours.
- j. Add a second aliquot of 3 µl of *DpnII* enzyme at the end of the day and shake overnight.
- k. Add a further 3 µl of *DpnII* enzyme the next morning and shake for at least three more hours.

3. Ligation (1 day)

- a. Place the digests on the 65 °C block for 20 min to heat inactivate the restriction enzyme.
- b. Cool the digest on ice and cool the ThermoMixer to 16 °C.
- c. Make Ligation Mix (make a master mix for multiple libraries):

PCR grade water	125 µl
10x ligation buffer	36 µl
- d. Add 161 µl Ligation Mix to digest.
- e. Add 4 µl T4 ligase and shake at 300 rpm/16 °C for ~22 h.

4. De-crosslinking (1 day)

Add 2 µl Proteinase K (20 mg/ml) and incubate at 65 °C overnight.

5. DNA extraction (1 day)

- a. Add 4 µl of RNase to the ligation reaction and incubate at 37 °C for 30 min.
- b. Add 350 µl phenol-chloroform-isoamyl alcohol and vortex thoroughly.
- c. Transfer to a 2 ml PLG Light tube and spin for 5 min/15,000 x g/RT.
- d. Transfer the upper layer, ~350 µl, to a new tube.
- e. Add 35 µl of 3 M NaOAc and 2 µl of GlycoBlue and mix by inversion.
- f. Add 960 µl of 100% ethanol and mix by inversion.
- g. Freeze at -20 °C overnight.

SAFE STOPPING POINT: DNA should be precipitated for at least one night; this can be extended to several days without affecting recovery.

- h. Pre-cool centrifuge to 4 °C during incubation.
- i. Centrifuge for 45 min/21,000 x g/4 °C and carefully discard liquid.
- j. Wash pellet with 1 ml 70% ice-cold ethanol.
- k. Spin for 10 min/21,000 x g/4 °C.
- l. Remove ethanol and repeat ethanol wash (steps A5j-A5k).
- m. Remove ethanol, spin in a microfuge and use a pipette to remove residual ethanol.
Dry at room temperature and resuspend pellet. When proceeding with LI-Capture-C, resuspend in 124 µl water. For Tag-Capture-C, the amount of water depends on the input, as several tagmentation reactions are set up in parallel (explained below). Use 20 µl water per multitude of 15,000 cells and round up (e.g., 30,000 cells:40 µl water; 50,000 cells:80 µl water).

SAFE STOPPING POINT: Store at -20 °C.

6. Quantification and QC (3 h)

Note: The NanoDrop or similar spectrophotometers are unreliable for quantifying 3C libraries due to residual DTT from the ligation buffer.

- a. Assess 1 µl of 3C library using the Agilent Genomic ScreenTape system to determine the size distribution of library ([Appendix I](#)).
- b. For LI-Capture-C—take 3 µl of 3C library and 27 µl water to make a 1/10 dilution; for Tag-Capture-C—take 1 µl of 3C library and 29 µl water to make a 1/30 dilution.

Perform qPCR for digestion efficiency, using the dilution. Genomic DNA can be used as the undigested control ([Appendix I](#)). Digestion efficiency should be > 75%.

B. Addition of Illumina sequencing adaptors by ligation or tagmentation to generate Capture-C libraries

LI-Capture-C

–Based on the protocol for the NEBNext Ultra DNA Library Prep Kit for Illumina (E7370)

Notes:

- a. *In this step, it is important to maintain library complexity by maximizing input DNA and minimizing DNA losses during the reactions and cleanups. When more than 1 µg of DNA is available, perform parallel preparations. These can be pooled during cleanups and at the end.*
- b. *The NEB protocol often recommends removing the DNA from the beads with a few µl more water than is necessary for the next step in the protocol. One can avoid doing this to minimize losses, but this means being very careful with the bead cleanups as contamination by the beads or ethanol can inhibit the following reaction.*
- c. *We use the Covaris S220 Focused Ultrasonicator for sonication. If using a different model sonicator, use high molecular weight gDNA to optimize sonication for a modal distribution around 200 bp in size ([Appendix II](#)).*
- d. *The Illumina indices that you put on to the library need to be complementary to the Nimblegen HE blocking oligonucleotides required in the capture step (the kits do not match; [Appendix II](#)). Buying both primer set 1 and primer set 2 eliminates waste of the more expensive Nimblegen blocking kit.*
- e. *No size selection is necessary as adaptor dimers will not be captured. All Ampure XP bead cleanups are performed with 1.8x volumes of beads.*

1. Sonication (1 h)

- a. Take 120 µl of 3C library (up to 1 µg) and transfer to a Covaris microtube. Avoid making bubbles.
- b. Shear DNA to 200 bp with the following settings:
 - Duty cycle: 10%
 - Intensity: 5
 - Cycles per burst: 200
 - Time: 360 sec
 - Set mode: Frequency sweeping
- c. Ampure XP SPRI bead cleanup:
 - i. Transfer reaction from the Covaris microtube to a lowBind DNA 1.5 ml microtube.
 - ii. Add 220 µl beads, pipette up and down 10 times, allow to bind at RT for 5 min.
 - iii. Place on the DynaMag, discard liquid when clear.

- iv. Add 500 μ l of fresh 80% ethanol without removing from the DynaMag. Avoid disturbing beads by running the ethanol down the front of the tube. Allow to sit on the rack for 30 sec at RT, then remove ethanol.
 - v. Remove ethanol and add another 500 μ l of fresh 80% ethanol, allow to sit for 30 sec at RT and remove all the ethanol.
 - vi. Spin down on microfuge and remove residual ethanol with a P10 pipette.
 - vii. Air dry at room temperature on the DynaMag until matt in appearance.
Note: Be careful not to over dry the beads as this will result in increased DNA losses.
 - viii. Resuspend beads in 56 μ l water to elute DNA, mix by pipetting 10 times.
 - ix. Incubate at RT for 2 min.
 - x. Recover 55.5 μ l for End Repair reaction.
- d. Optional: Assess 1 μ l of sonicated material using the Agilent D1000 (Regular or High Sensitivity, depending on your input) ScreenTape system ([Appendix II](#)).
- SAFE STOPPING POINT: Store at -20 °C.*
2. End repair and adaptor ligation (2 h)
- a. Combine 55.5 μ l of DNA (5 ng-1 μ g), 6.5 μ l 10x End Repair Reaction Buffer (Green) and 3 μ l End Prep Enzyme Mix (Green). Mix by pipetting and briefly spin down.
 - b. Incubate in a thermocycler:
20 °C for 30 min
65 °C for 30 min
 - c. Add 15 μ l Blunt/TA Ligase Master Mix (Red), 2.5 μ l NEBNext Adaptor* (Red), 1 μ l Ligation Enhancer (Red). Mix by pipetting and briefly spin down.
**Note: If DNA input is < 100 ng, dilute the NEBNext Adaptor for Illumina (provided at 15 μ M) 10-fold in 10 mM Tris-HCl or 10 mM Tris-HCl with 10 mM NaCl to a final concentration of 1.5 μ M, use immediately.*
 - d. Incubate for 15 min at 20 °C in a thermocycler.
 - e. Add 3 μ l of USER Enzyme (Red). Mix by pipetting.
 - f. Incubate for 15 min at 37 °C in a thermocycler.
 - g. Ampure XP SPRI bead cleanup:
 - i. Transfer reaction to a lowBind DNA 1.5 ml microtube.
 - ii. Add 160 μ l beads, pipette up and down 10 times, allow to bind at RT for 5 min.
 - iii. Place on the DynaMag, discard liquid when clear.
 - iv. Add 500 μ l of fresh 80% ethanol without removing from the DynaMag. Avoid disturbing beads by running the ethanol down the front of the tube. Allow to sit on the rack for 30 sec at RT, then remove ethanol.
 - v. Remove ethanol and add another 500 μ l of fresh 80% ethanol, allow to sit for 30 sec at RT and remove all the ethanol.
 - vi. Spin down in a microfuge and remove residual ethanol with a P10 pipette.
 - vii. Air dry at room temperature on the DynaMag until matt in appearance.

Note: Be careful not to over dry the beads as this will result in increased DNA losses.

- viii. Resuspend beads in 29 μ l water, mix by pipetting 10 times.
- ix. Incubate at RT for 2 min.
- x. Transfer 28.5 μ l eluted DNA to a PCR tube.

3. PCR addition of indices (2 h)

Note: Carefully choose your index to allow for multiplexing using the Illumina pooling guidelines.

a. Perform PCR as described below.

Adaptor ligated 3C library	28.5 μ l	Step 1	98 °C	30 sec
NEB Universal primer	5 μ l	Step 2	98 °C	10 sec
NEB Index primer	5 μ l	Step 3	65 °C	30 sec
Herculase II Reaction Buffer (5x)	10 μ l	Step 4	72 °C	30 sec
dNTP Mix (25 mM)	0.5 μ l	Step 5	Go to Step 2	9 cycles
Herculase II Fusion DNA Polymerase	1 μ l	Step 6	72 °C	5 min
	Total 50 μ l	Step 7	4 °C	Hold

b. Ampure XP SPRI bead cleanup:

- i. Transfer reaction to a lowBind DNA 1.5 ml microtube.
- ii. Add 90 μ l beads, pipette up and down 10 times, allow to bind at RT for 5 min.
- iii. Place on the DynaMag, discard liquid when clear.
- iv. Add 500 μ l of fresh 80% ethanol without removing from the DynaMag. Avoid disturbing beads by running the ethanol down the front of the tube. Allow to sit on the rack for 30 sec at RT, then remove ethanol.
- v. Remove ethanol and add another 500 μ l of fresh 80% ethanol, allow to sit for 30 sec at RT and remove all the ethanol.
- vi. Spin down in a microfuge and remove residual ethanol with a P10 pipette.
- vii. Air dry at room temperature on the DynaMag until matt in appearance.
Note: Be careful not to over dry the beads as this will result in increased DNA losses.
- viii. Resuspend beads in 31 μ l water, mix by pipetting 10 times.
- ix. Incubate at RT for 2 min.
- x. Recover 30 μ l.

c. Assess 1 μ l of material using the Agilent D1000 ScreenTape system ([Appendix II](#) page 3).

d. Quantify library using Qubit dsDNA BR assay kit.

SAFE STOPPING POINT: Store at -20 °C.

Tag-Capture-C

Notes:

- a. In this step, it is important to maintain library complexity by maximizing input DNA and minimizing DNA losses during the reactions and cleanups. The maximum input for a tagmentation reaction is 50 ng. When input is higher, run several separate reactions in parallel; these can be pooled after the cleanup.

- b. We use custom designed index primers, for which the sequences can be found in the Appendix. We order these as TruGrade Primers from IDT.
- c. No size selection is necessary as adaptor dimers won't be captured. All Ampure XP bead cleanups are performed with 1.8x volumes of beads.

1. Tagmentation (1 day)

- a. Prepare tagmentation reactions with reagents from the Nextera DNA Sample Preparation Kit in a PCR tube:
 - < 50 ng 3C library in 20 µl water
 - 25 µl TD Buffer
 - 5 µl of TD Enzyme
- b. Mix by pipetting, followed by a quick spin.
- c. Incubate at 55 °C for 12-24 h to ensure all material is tagmented.
- d. Zymo spin column cleanup:
 - i. Label tubes and add 180 µl of Zymo DNA Binding Buffer.
 - ii. Add 50 µl of tagmented DNA. Gently pipette up and down 10 times to mix.
 - iii. Transfer sample mixture to a Zymo spin column in a collection tube.
 - iv. Centrifuge for 30 sec at 15,000 x g. Discard the flow-through.
 - v. Add 200 µl Zymo DNA Wash Buffer to the column. Centrifuge for 30 sec at 15,000 x g.
 - vi. Repeat step B1d.v (Tag-Capture-C section) and centrifuge in addition 60 sec at 15,000 x g.
 - vii. Transfer the column to a 1.5 ml tube and add 34.5 µl of water to the column matrix. Incubate at RT for 2 min.
 - viii. Centrifuge for 1 min at 15,000 x g to elute the DNA.
- e. Assess 1 µl of material using the Agilent High Sensitivity D1000 ScreenTape system ([Appendix III](#)).

2. PCR addition of indices (2 h)

Note: Carefully choose your index to allow for multiplexing using the Illumina pooling guidelines.

- a. Perform PCR as described below.

Tagmented 3C library	33.5 µl
Universal primer (25 µM)	2 µl
Index primer (25 µM)	2 µl
KAPA HiFi Fidelity Buffer (5x)	10 µl
KAPA dNTP Mix (10 mM)	1.5 µl
KAPA HiFi DNA Polymerase (1 U/µl)	1 µl
	Total 50 µl

Step 1	72 °C	5 min
Step 2	98 °C	30 sec
Step 3	98 °C	10 sec
Step 4	63 °C	30 sec
Step 5	72 °C	1 min
Step 6	Go to Step 3	9 cycles
Step 7	72 °C	3 min
Step 8	4 °C	Hold

- b. Ampure XP SPRI bead cleanup:
 - i. Transfer reaction to a lowBind DNA 1.5 ml microtube.
 - ii. Add 90 µl beads, pipette up and down 10 times, allow to bind at RT for 5 min.

- iii. Place on the DynaMag, discard liquid when clear.
 - iv. Add 500 μ l of fresh 80% ethanol without removing from the DynaMag. Avoid disturbing beads by running the ethanol down the front of the tube. Allow to sit on the rack for 30 sec at RT, then remove ethanol.
 - v. Remove ethanol and add another 500 μ l of fresh 80% ethanol, allow to sit for 30 sec at RT and remove all the ethanol.
 - vi. Spin down in a microfuge and remove residual ethanol with a P10 pipette.
 - vii. Air dry at room temperature on the DynaMag until matt in appearance.
Note: Be careful not to over dry the beads as this will result in increased DNA losses.
 - viii. Resuspend beads in 31 μ l water, mix by pipetting 10 times.
 - ix. Incubate at RT for 2 min.
 - x. Recover 30 μ l.
- c. Assess 1 μ l of material using the Agilent D1000 ScreenTape system ([Appendix III](#) page 4).
 - d. Quantify library using Qubit dsDNA BR assay kit.
- SAFE STOPPING POINT: Store at -20 °C.**

C. Enrichment for viewpoints of interest by oligonucleotide capture

–Based on the Nimblegen SeqCap SR User’s Guide Chapter 5-7

Note: The blocking oligonucleotides must be complementary to the adapters in your library. If performing LI-Capture-C, you can use the commercially available blocking oligonucleotides from Nimblegen. If you are following the Tag-Capture-C protocol, you can contact IDT and design blocking oligonucleotides using their xGen Custom Blocking Oligos service. We used a universal blocker for both the i5 adapter and i7 adapters, so we only needed 2 blocking oligonucleotides for all possible index combinations.

1. Hybridization preparation (1 day)

Note: Ensure capture oligonucleotides for different experiments are kept apart. The oligonucleotides are hugely in excess and tiny amounts of contamination can lead to spurious results. We would recommend that you do not order oligonucleotides that should not be mixed from the manufacturer at the same time. Avoid having buffers and blocking oligonucleotides out at the same time as the capture oligonucleotides. Oligonucleotides that are going to be pooled for a single experiment, can be ordered together. Consider ordering large designs in as a pooled set (at equimolar concentrations).

Biotinylated capture oligonucleotides

- a. Reconstitute individual or pools of oligonucleotides to a stock concentration so that each oligo is stored at 2.9 μ M (10 nmol in 3.46 μ l of PCR grade water)–or any high concentration.
- b. Generate pools of oligonucleotides by mixing in exact 1:1 stoichiometric ratio.

Nimblegen blocking oligonucleotides (LI-Capture-C)

Resuspend Nimblegen oligonucleotides according to the Nimblegen protocol:

- a. Spin down tubes in a minifuge.

- b. Add 120 μ l PCR grade water to the HE Universal Oligo tube (1 mM); vortex and spin briefly.
- c. Add 10 μ l PCR grade water to the HE Index Oligo tube (1 mM); vortex and spin briefly.
- d. Add 480 μ l PCR grade water to the Post-LM-PCR Oligos; vortex and spin briefly.

Custom blocking oligonucleotides (Tag-Capture-C)

Resuspend the blocking oligonucleotides according to manufacturer's instructions.

Multiplexing of samples

Mix up to 1 μ g of each differentially indexed sample at exactly 1:1 ratios by mass to generate a multiplexed library.

2. Hybridization (3 days)

- a. Heat vacuum centrifuge to 50 °C.
- b. Prepare hybridization reaction for the number of pooled libraries—up to 6 libraries are captured in a single tube. For more libraries, a master mix may be made in one tube and divided into multiple tubes.

5 μ g COT DNA (5 μ l of stock)/library

1 nmol TS-HE Universal Oligo (1 μ l of 1 mM stock aliquot)/library

1 nmol of TS-HE Index Oligos (1 μ l of 1 mM stock for each index)

up to 1 μ g of each uniquely indexed 3C library

Note: For fewer than 100,000 cells you may not recover > 1 μ g after indexing. Use as much as possible.

- c. Vacuum centrifuge at 50 °C with tube lids open (rather than pierced) until sample is completely dry. Avoid drying for a long time after the liquid is gone.
- d. For each library in the capture reaction add:
 - 7.5 μ l 2x hybridization buffer (vial 5)
 - 3 μ l hybridization component A (vial 6)
- e. Carefully reconstitute the DNA by pipetting and vortexing followed by briefly spinning down.
- f. **Replace buffers and blocking oligonucleotides in freezer prior to proceeding.**
- g. Preheat the Eppendorf ThermoMixer to 95 °C.
- h. Preheat a PCR block to 47 °C, the lid should be heated to 57 °C.
- i. Heat 4.5 μ l per library of pooled biotinylated oligonucleotide probes to 47 °C in a PCR tube.
Note: The capture will be at 47 °C for ~72 h. A high-quality PCR tube is required to avoid sample loss through evaporation.
- j. Denature the 3C library and blocking oligonucleotides by heating to 95 °C for 10 min.
- k. After 10 min, quickly spin the denatured library and add the entire volume to the pooled biotinylated oligonucleotides without removing from the thermocycler.
- l. Mix and spin briefly before replacing in the 47 °C PCR machine for 64-72 h.

3. Washing and recovery of captured material (1 day)

Notes:

- i. *Heat reagents on an Eppendorf ThermoMixer as it is more reliable than a water bath and allows the samples to be shaken.*

- ii. For capture of multiplexed libraries, scale the beads and wash buffers accordingly to the number of libraries.
 - iii. The Stringent Wash buffer, Wash Buffer I and Bead Wash Buffer are quite voluminous and may need to be made in multiple tubes. Split Wash Buffer I 1:2 as one third is heated.
- a. Dilute the wash buffers as per table below which allows a slight excess of each.
 - i. Place the Stringent Wash Buffer at 47 °C on the ThermoMixer.
 - ii. Place 100 µl Wash Buffer I per captured library with an additional 20 µl at 47 °C on the ThermoMixer.

Number of libraries in capture:	Buffer to add for 1x			Water to add for 1x		
	1	2	3	1	2	3
Stringent Wash Buffer (vial 4)	42	82	122	378	738	1,098
Wash Buffer I (vial 1)	34	64	94	306	576	846
Wash Buffer II (vial 2)	22	42	62	198	378	558
Wash Buffer III (vial 3)	22	42	62	198	378	558
Bead Wash Buffer (vial 7)	210	410	610	315	615	915

Note: The streptavidin beads will likely stick to the interior of the tube. This is particularly a problem with some makes of low bind tubes. Use tubes with minimal adhesion. We use Safeseal tubes (Sorenson) for this step.

- b. Prepare the streptavidin beads (M270):
 - i. Allow the beads to heat to RT for 30 min prior to use.
 - ii. Aliquot 100 µl per captured library into a single 1.5 ml Eppendorf.
 - iii. Place on the DynaMag and remove liquid once clear.
 - iv. Add 200 µl of 1x Bead Wash Buffer per captured library and vortex to resuspend the beads, spin briefly.
 - v. Replace on the DynaMag for 5 min and remove liquid once clear.
 - vi. Repeat wash steps (C3b.iii-C3b.v) for a total of two washes.
 - vii. Resuspend the beads in their original volume (*i.e.*, 100 µl per captured library) with Bead Wash Buffer (1x) and aliquot into an appropriate number of tubes 1.5 ml Eppendorf tubes.
 - viii. Place on the DynaMag, remove and discard the liquid only when ready to add the capture sample from the thermocycler. **Do not allow beads to dry out.**
- c. Binding of biotinylated oligonucleotides:
 - i. Transfer the hybridization reactions to the streptavidin beads and mix thoroughly by pipetting 10 times.
 - ii. Spin briefly if necessary to pool all of the sample in the bottom of the tube.
 - iii. Place on the ThermoMixer at 47 °C, 600 rpm for 45 min. Mix the beads with a pipette every 15 min to avoid beads settling at the bottom of the tube.

Note: Steps C3d.i-C3d.vii should be done quickly to maintain the temperature at 47 °C.

- d. Washing the streptavidin beads and bound DNA:
- i. Add 100 μ l of heated Wash Buffer I (47 °C) per captured library to the beads and bound DNA.
 - ii. Mix by vortexing.
 - iii. Place on the DynaMag, carefully discard all the supernatant when clear.
 - iv. Add 200 μ l Stringent Wash buffer per captured library heated to 47 °C and mix.
 - v. Incubate for 5 min at 47 °C.
 - vi. Place the tubes on the DynaMag, carefully discard the supernatant when clear.
 - vii. Repeat stringent wash (C3d.iv-C3d.vi) for a total of 2 washes.
 - viii. Add 200 μ l Wash Buffer I (1x) at RT per captured library and mix by vortexing for 2 min. Spin briefly to ensure no sample is lost in the lid.
 - ix. Place the tubes on the DynaMag, carefully discard the supernatant when clear.
 - x. Add 200 μ l Wash Buffer II (1x) per captured library and mix by vortexing for 1 min. Spin briefly to ensure no sample is lost in the lid.
 - xi. Place the tubes on the DynaMag, carefully discard the supernatant when clear.
 - xii. Add 200 μ l Wash Buffer III (1x) and mix by vortexing for 30 sec. Spin briefly to ensure no sample is lost in the lid.
 - xiii. Place the tubes on the DynaMag, carefully discard the supernatant when clear.
 - xiv. Remove from the DynaMag and resuspend beads in 40 μ l PCR grade water per capture library. Do not discard beads–DNA is not eluted but amplified off the beads.
 - xv. Store at -20 °C or proceed to PCR amplification.
- SAFE STOPPING POINT: Store at -20 °C.*

4. PCR amplification of captured material (1 day)

- a. Amplify the captured fragments in two separate reactions per captured library in a total of 14 cycles.

Hybridized beads	20 μ l
KAPA HiFi HotStart ReadyMix	25 μ l
POST-LM_PCR Oligo 1&2	5 μ l
	Total 50 μ l

Step 1	98 °C	45 sec
Step 2	98 °C	15 sec
Step 3	60 °C	30 sec
Step 4	72 °C	30 sec
Step 5	Go to Step 2	14 cycles
Step 6	72 °C	60 sec
Step 7	4 °C	Hold

Note: Only use 20 μ l of DNA/Bead–using 40 μ l has the potential to inhibit the reaction. Two separate PCR reactions can be performed or 20 μ l can be stored at -20 °C.

- b. Ampure XP SPRI bead cleanup:
- i. Transfer reactions to lowBind DNA 1.5 ml microtubes. Combine up to six reactions.
 - ii. Add 90 μ l beads per reaction, pipette 10 times to mix, allow to bind at RT for 5 min.
 - iii. Place on the DynaMag, discard liquid when clear.

- iv. Add 500 μ l of fresh 80% ethanol without removing from the DynaMag. Avoid disturbing beads by running the ethanol down the front of the tube. Allow to sit on the rack for 30 sec at RT, then remove ethanol.
 - v. Remove ethanol and add another 500 μ l of fresh 80% ethanol, allow to sit for 30 sec at RT and remove all the ethanol.
 - vi. Spin down in a microfuge and remove residual ethanol with a P10 pipette.
 - vii. Air dry at room temperature on the DynaMag until matt in appearance.
Note: Be careful not to over dry the beads as this will result in increased DNA losses.
 - viii. Resuspend beads in 51 μ l water, mix by pipetting 10 times.
 - ix. Incubate at RT for 2 min.
 - x. Recover 50 μ l.
- c. Assess 1 μ l of material using the Agilent D1000 ScreenTape system to confirm same profile as input material.
 - d. Quantify with Qubit BR Kit.
SAFE STOPPING POINT: Store at -20 °C.
5. Double capture (3 days)
Note: Capture efficiency is improved by 100-1,000 fold by performing a second capture step. This increases the proportion of captured reads from 1% to 50%. This means only ~1 million reads are needed per viewpoint per sample and thus reduces sequencing requirements. Therefore, double capture is standard practice in our lab.
 - a. Pool amplified material of the parallel PCR reactions.
 - b. Prepare the second capture hybridization reaction using all captured material (up to 2 μ g) in a single reaction. Combine:
 - i. 5 μ g COT DNA (5 μ l of stock)
 - ii. 1 nmol TS-HE Universal Oligo (1 μ l of 1 mM stock aliquot)
 - iii. 1 nmol of TS-HE Index Oligos (1 μ l of 1 mM stock in total; so if you have 4 samples, use 0.25 μ l of each index for a total of 1 μ l; consider making a dilution if you have many samples)
 - iv. Up to 2 μ g of amplified captured material
 - c. Perform hybridization, washes, and amplification as described with a 24-h hybridization at 47 °C and treating the material as a single library.
SAFE STOPPING POINT: Store at -20 °C.
 6. Library quantification and sequencing (1 day)
 - a. Use qPCR with the KAPA Library Quantification Kit to calculate the concentration of adaptor containing fragments.
Make 1:10,000 and 1:20,000 dilutions of the captured material for quantification.
 - b. Sequence using Illumina paired-end sequencing with either 150 bp or 75 bp reads (300-cycle and 150-cycle kits).

Data analysis

The data analysis of low-input Capture-C experiments has been described (Oudelaar *et al.*, 2017). Scripts are available in the GitHub repository (<https://github.com/oudelaar/CaptureC/>).

Recipes

1. Fresh lysis buffer, cool to 4 °C—for 10 samples (5 ml):
 - 10 mM Tris pH 8 (50 µl 1 M)
 - 10 mM NaCl (10 µl of 5 M)
 - 0.2% Igepal CA-630 (100 µl of 10%)
 - 1x cComplete Protease Inhibitor Cocktail (200 µl of 25x)
 - 4.64 ml PCR grade water

Acknowledgments

The low-input Capture-C protocols are optimized adaptations from the Next-Generation Capture-C method (Davies *et al.*, 2016). We thank the Wellcome Trust (Wellcome Trust Genomic Medicine and Statistics PhD Programme, reference 105281/Z/14/Z; Wellcome Trust Strategic Award, reference 106130/Z/14/Z) and the Medical Research Council (MRC Core Funding and Centenary Award reference 4050189188) for funding our work. The authors declare that there is no conflict of interest.

References

1. Buenrostro, J. D., Giresi, P. G., Zaba, L. C., Chang, H. Y. and Greenleaf, W. J. (2013). [Transposition of native chromatin for fast and sensitive epigenomic profiling of open chromatin, DNA-binding proteins and nucleosome position.](#) *Nat Methods* 10(12): 1213-1218.
2. Davies, J. O., Oudelaar, A. M., Higgs, D. R. and Hughes, J. R. (2017). [How best to identify chromosomal interactions: a comparison of approaches.](#) *Nat Meth* 14: 125-134.
3. Davies, J. O., Telenius, J. M., McGowan, S. J., Roberts, N. A., Taylor, S., Higgs, D. R. and Hughes, J. R. (2016). [Multiplexed analysis of chromosome conformation at vastly improved sensitivity.](#) *Nat Methods* 13(1): 74-80.
4. Dekker, J., Rippe, K., Dekker, M. and Kleckner, N. (2002). [Capturing chromosome conformation.](#) *Science* 295(5558): 1306-1311.
5. Oudelaar, A. M., Davies, J. O. J., Downes, D. J., Higgs, D. R. and Hughes, J. R. (2017). [Robust detection of chromosomal interactions from small numbers of cells using low-input Capture-C.](#) *Nucleic Acids Res* 22(45): e184.